

vista®

for...

TRANSPARENT PLASTIC SHEET




Vista® simultaneously measure transmittance haze% and visible-range transmittance color of transparent plastic sheet in a single measurement. Vista® includes easy to navigate embedded EasyMatch® Essentials software, providing all the scales and indices important in the plastics industry including CIE L*a*b*, YI E313, WI 313, and haze% in conformance with ASTM D1003- Procedure B.

With modern data communications and a capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of transparent plastic sheet samples.



Simultaneously measure the transmission color and haze% of transparent sheet plastic

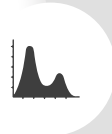

Full spectrum, balanced LED illumination for years of worry-free operation




Operator friendly one-touch standardization – no accessories required



Embedded EasyMatch® Essentials software – no PC required



Smart user interface and data management

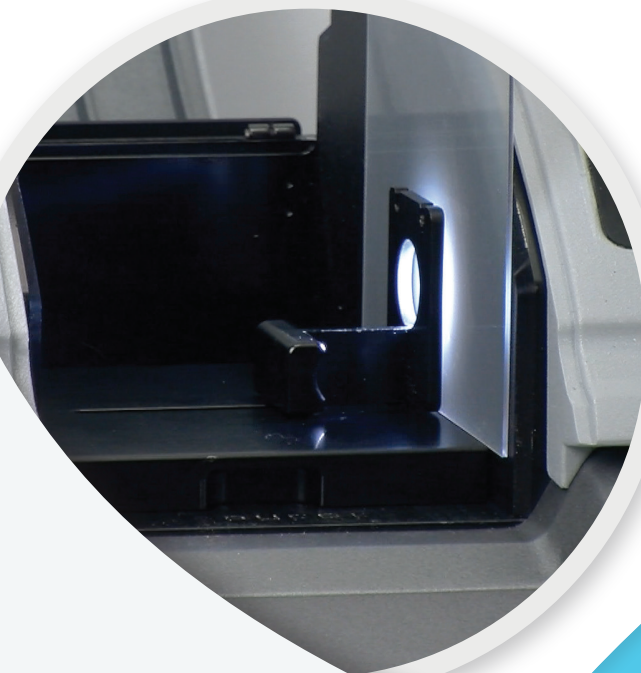



measure
color
the way your
eye sees it



The Best Solution for Color Measurement of Transparent Plastic Sheet


Simultaneously measures the transmittance color and haze of transparent preforms

Get accurate haze and visible-range color results in a single measurement. Haze conforms to ASTM D1003 Procedure B.

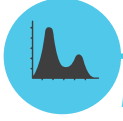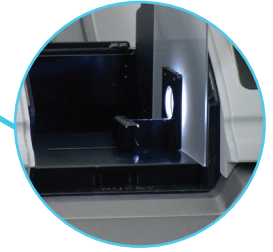


Full spectrum, balanced LED illumination

Illuminating from 400-700nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements.

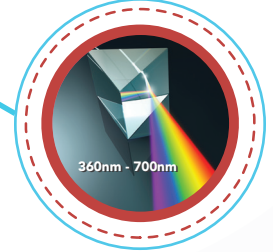

One-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

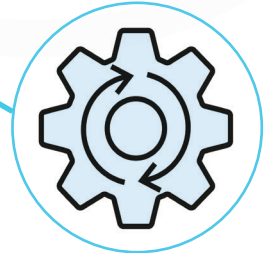

Embedded EasyMatch® Essentials software

The onboard EasyMatch® Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard transmittance scales and indices widely used in industry.



Smart user interface and data management - no PC required

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.

Keys for the successful measurement of Transparent Plastic Sheet

The following are required to measure Transparent Plastic Sheet with Vista®:

ITEM	HunterLab PN#
Vista® Transmittance Spectrophotometer	Vista®
Multi-Function Sample Holder	D02-1017-192

Multi-Function Sample Holder

