


**Compact
Powerful
Precise**

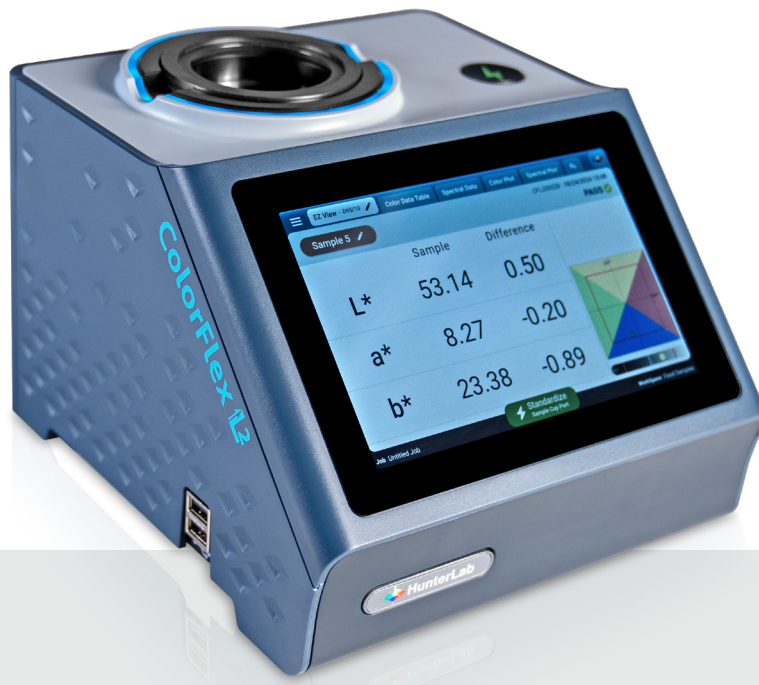

ColorFlex[®] L2

*Level UP your Color Confidence with the new
ColorFlex L2 Spectrophotometer.*

As simple as you want, as powerful as you need.

Elevate your color accuracy and user experience like never before. Designed to streamline and expedite your quality control processes, ColorFlex L2 Spectrophotometer gives you unsurpassed confidence and trust in your outcomes. Compact and ergonomic build requires minimal space while offering enhanced durability and reliability.


Explore the Next-Gen ColorFlex®...L2



Precise color measurement for everything from opaque solids, liquids, powders, granules and pellets to translucent solids and liquids. We have upgraded the most relied upon 45°/0° spectrophotometer in the world, ColorFlex EZ, to the new ColorFlex L2 with the most advanced color technology.



Explore all the additional features and benefits below.

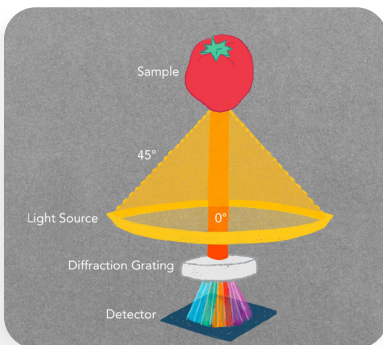

45°a:0° Annular Illuminated Optical Geometry

Superior Consistency and Color Accuracy: Enhancing the efficiency of Quality Control workflows and delivery of superior quality products repeatedly.



Compact with Sealed Spillproof Case

Space-Saving and Easy Maintenance: Occupies minimal space and is designed to withstand the demands of busy lab environments.

Powerful Standalone Functionality

No computer required: With huge onboard processing and storage capacity, ColorFlex L2 can function as an independent color workstation. Connect monitor, keyboard, and mouse if desired.




Sliced Cheese SPICES flour POWDERED SPORTS DRINKS Yogurt
Calcium Carbonate Puree mustards barium sulfate
BABY FOOD Coffee - Green beans, Roasted & Ground CEMENT
Baking Powders Citrus Juice Herbs Cocoa Powder



Modern Tablet User Interface with Touchscreen & Color Display

Durable and Efficient: Facilitates accurate color measurement in seconds, with a robust industrial glass design suitable for harsh environments.

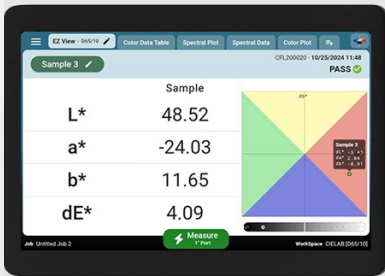


New & Improved Color Software

Essentials 2.0: Take precise and reliable color measurement in seconds with next-gen user-friendliness and seamlessly analyze multiple tristimulus and spectral data sets. It's that powerful.


Self-Guided and Intuitive Out-of-Box Startup

Time-Saving and User-Friendly: Ready to use straight out of the box with an integrated wizard for quick training, reducing operator time and easing concerns related to staff training and turnover.

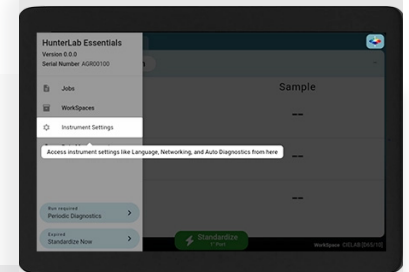

Built-in Camera for Sample Viewing and Screen Recording

Precision and Communication: Ensures accurate sample positioning to eliminate measurement discrepancies and allows easy sharing of recordings for swift quality control decisions.

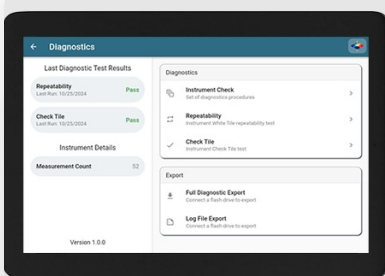

Built-in Advanced Diagnostics

On-Board Diagnostic Health Check: Ensures the HunterLab ColorFlex L2 is always operating at peak performance, continuing to build on the legendary history of reliability of the ColorFlex family of instruments.

Powdered Sugar	BAKING CHIPS	Concrete	Powdered Functional Beverages	Fertilizers
Magnesium Sulfate	tea drinks (powder)	Starch	shortening	WHITE SUGAR
RAW GRAINS	Brown Sugar	powdered cheese	SEASONINGS	Orange Juice



Advanced Data Export and Communication

Full suite of data communication options: Features USB ports, HDMI, and support for mouse and keyboard, enabling easy data export/printing.



Connectivity with LIMS and SPC systems

Simplifies QC workflows: Aids in trend analysis and workflow simplification across manufacturing sites.



CONTACT US

We've worked with thousands of companies, but you'll feel like our one and only. If color accuracy is crucial to your business, you've come to the right place—with the right people. Let's get in touch. Ask a question, Request a quote or just say hello!

OUR OFFICES

HunterLab

11491 Sunset Hills Road
Reston, VA 20190-5280
Website: www.hunterlab.com
Email: sales@hunterlab.com
Ph: 703-471-6870
Fax: 703-471-4237

Corporate Headquarters

HunterLab Europe

Christian Jansen, Managing Director
Dr.-August-Einsele-Ring 15
82418 Murnau, Germany
Website: www.hunterlab.de
Email: Christian.Jansen@hunterlab.com
Ph: +49 8841 9464

European Office

HunterLab Asia-Pacific

Leo Cheung, Regional Manager
Flat E, 13F Leahander Centre
28 Wang Wo Tsai Street Tsuen Wan
New Territories Hong Kong
Website: www.hunterlabcolor.cn
Email: sales-asiapacific@hunterlab.com
Ph: +852-2406 1723

Asia-Pacific Office

For more information, please contact HunterLab at 703-471-6870, sales@hunterlab.com or visit www.hunterlab.com